



# Product Specification Sheet

## NeoMag Slide Dock for EpiSep Slide™ (NSD I™)


**FOR RESEARCH USE ONLY:** This product is not for diagnostic or therapeutic use.

Doc.# ML081031-21

**CATALOG #:** A1101-1 (each)  
A1101-5 (5 Pack)


**DESCRIPTION:** The NSD I™ is a re-usable magnetic base for use with the EpiSep Slide® and EpiSep Multi-Purpose Slide®. The NSD I™ contains a single powerful neodymium permanent magnet that applies a focused magnetic field to the target cell capture area on the EpiSep Slide® or EpiSep Multi-Purpose Slide®. The EpiSep Slide® or EpiSep Multi-Purpose Slide® temporarily rests on the NeoMag Slide Dock during the addition of sample and reagents and is removed once the assay is complete for inspection or storage.

NSD I™

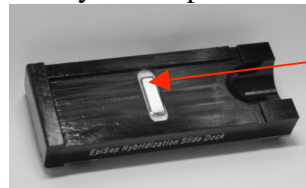

Neodymium Magnet

EpiSep Large Slide® mounted on the NSD I™

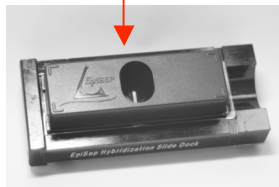

EpiSep Large Slide® is lifted off from the NSD I™



EpiSep Large Slide® is removed from the NSD I™



**SUPPLIED AS:** Package of 1 or 5 unit(s).

**STORAGE:** Store dry at room temperature.

**Designed for use with the EpiSep Standard Slide™, EpiSep Large Slide™ or the EpiSep Multi-Purpose Slide™ Only.**

### MATERIAL

**SAFETY DATA:** Contains plastic components. In the event of breakage some pieces may be sharp and can cause injury. Wear gloves and eye protection for clean up in the event of breakage.